



A. Synthesize tandem GAP-LOCK probes

1. Probe A
Target seq. A1 Linker seq. X
3'  5'
2. Probe B
Target seq. B1 Linker seq. Y
5'  3'

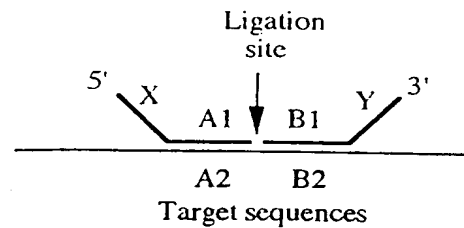
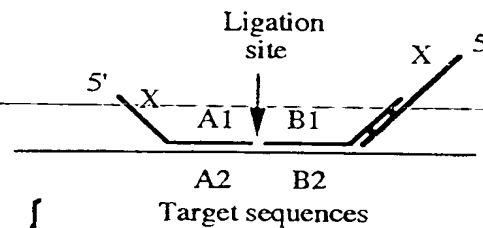
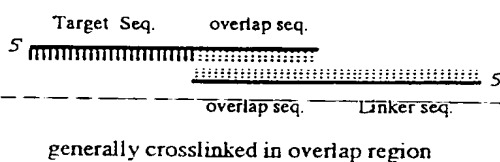
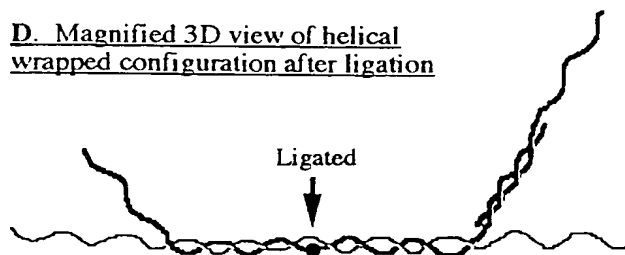
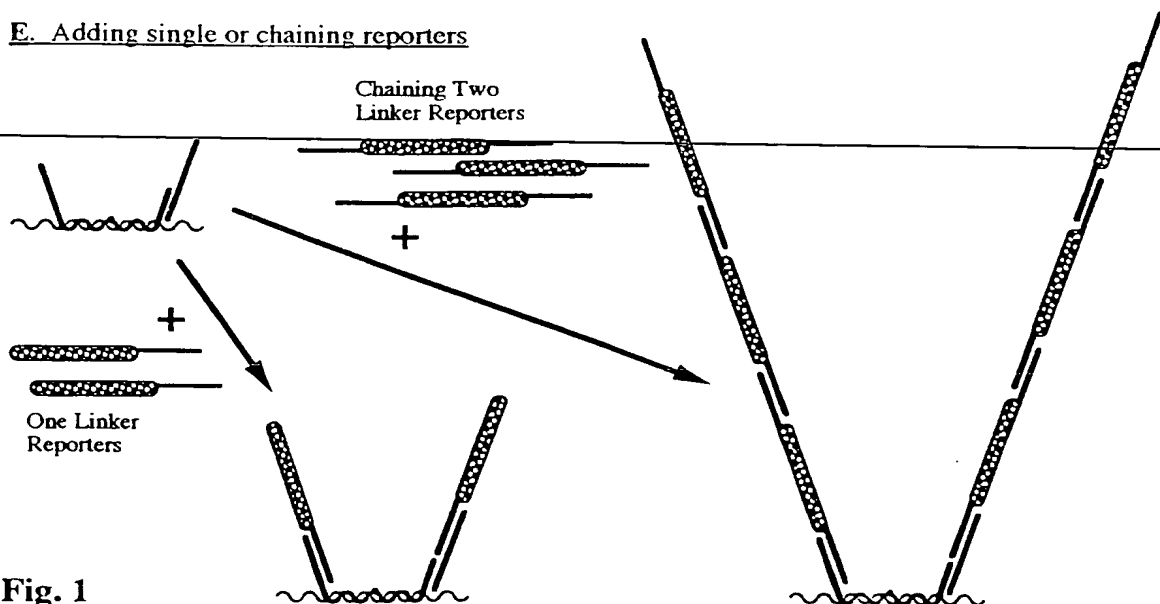
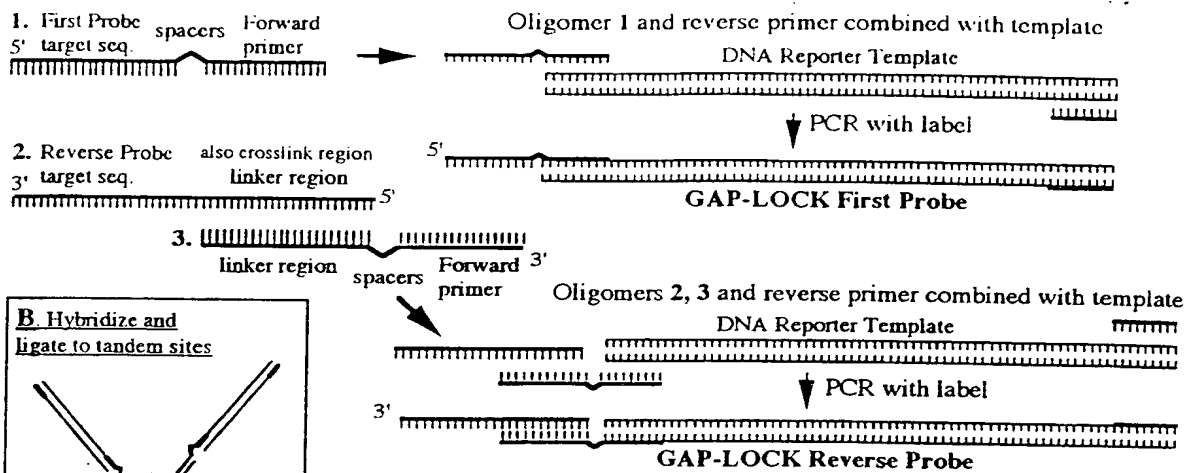
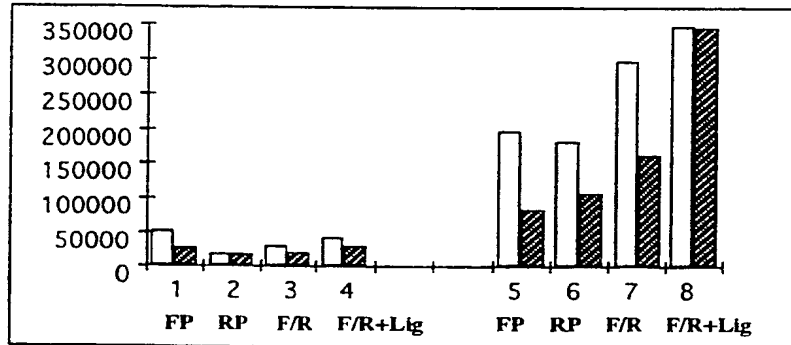
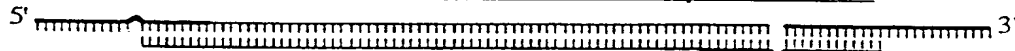
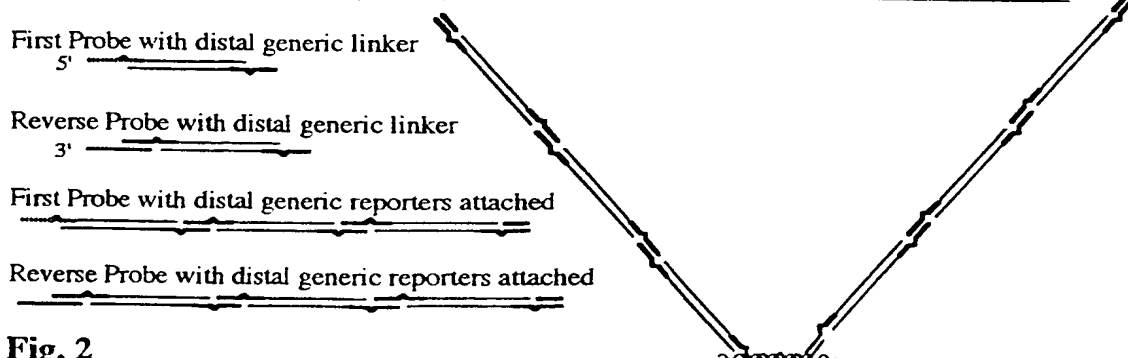
B. Hybridize Probes to target and ligate to form GAP-LOCKC. Alternatively, one probe made with reversing linker to provide same linker end as other probeD. Magnified 3D view of helical wrapped configuration after ligationE. Adding single or chaining reporters


Fig. 1

A. Alternate GAP-LOCK Method: with probes joined to reporters**C. Dot Blots of GAP-LOCK Probes: (Example 1)**



Bar 1-4 = 300 bp tails
 Bar 5-8 = 800 bp tails
 Open Bars = pre NaOH
 Stripe Bars = post NaOH
 F = First Probe
 R = Reverse Probe
 F/R = both probes
 F/R+Lig. = ligated both probes (First and Reverse)


D. Alternate GAP-LOCK with First Probe, Reverse Probe and Reporter combined**E. GAP-LOCK with generic reporters joined to distal linkers on the First and Reverse Probes****Fig. 2**

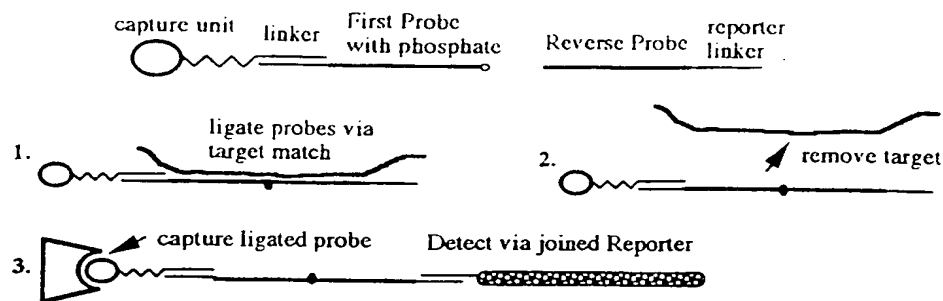Capture GAP-LOCK Method




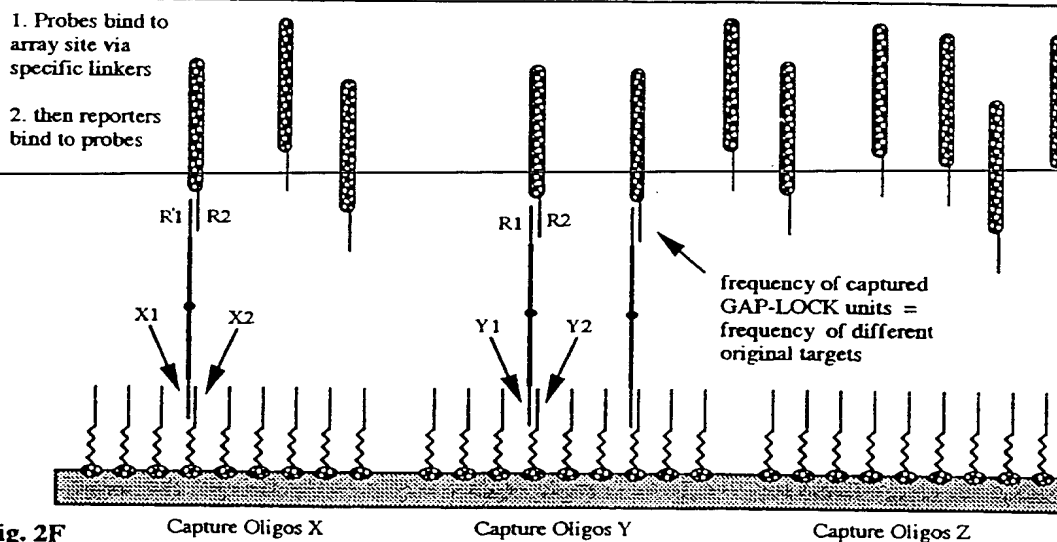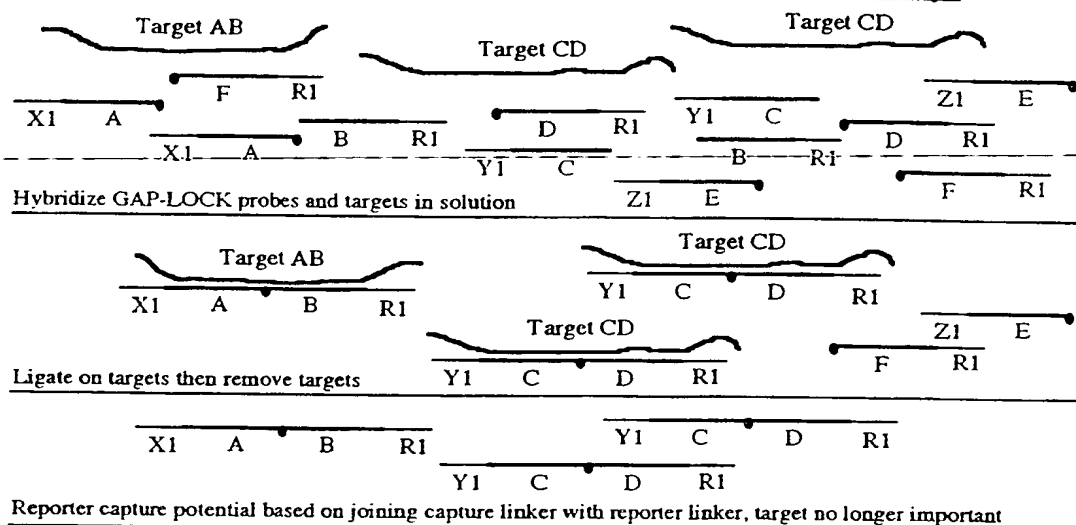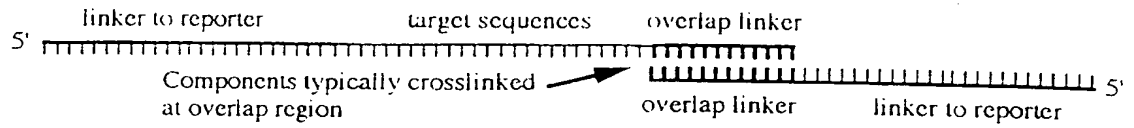Comparative detection of two or more targets with chips or microarrays





A. WRAP-PROBE Method: Component Assembly: Preferred EmbodimentB. Alternate WRAP-PROBE Embodiments

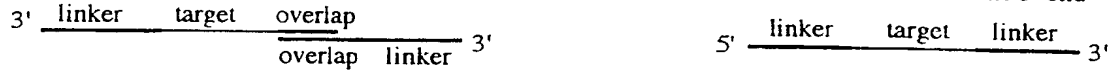Synthetic WRAP-PROBE made with 3' ends instead of 5' ends, or with one 5' end and one 3' end



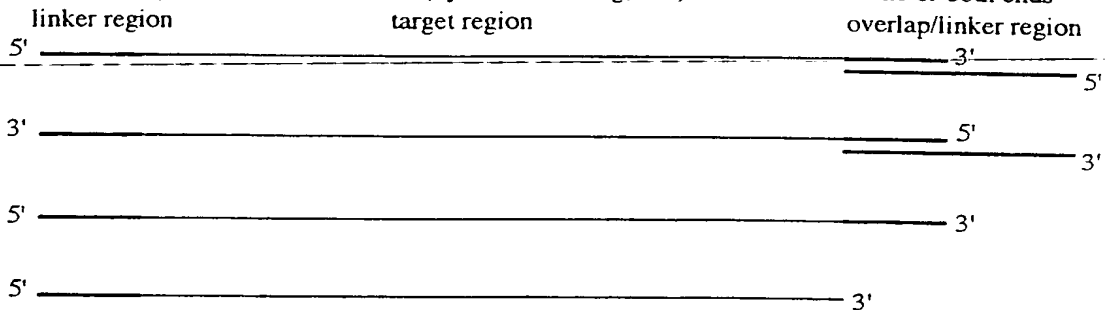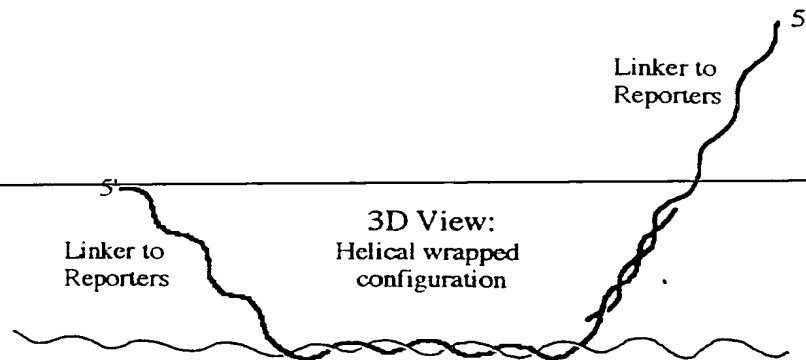
Enzymatically made WRAP-PROBE (by PCR, Cloning, etc.) with linkers on one or both ends

C. WRAP-PROBE Hybridized to Target Strand


The probe wraps around the target approximately one turn per each 10 bases of sequence.

Fig. 3

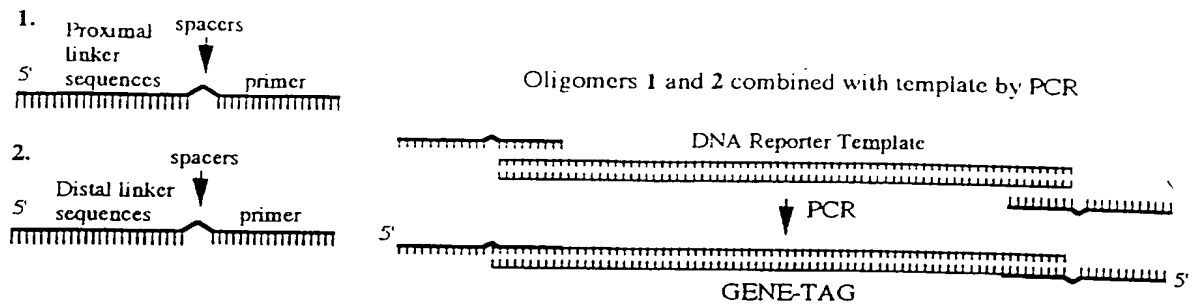A. Manufacture GENE-TAGs by oligomer synthesis and PCR with label or hapten



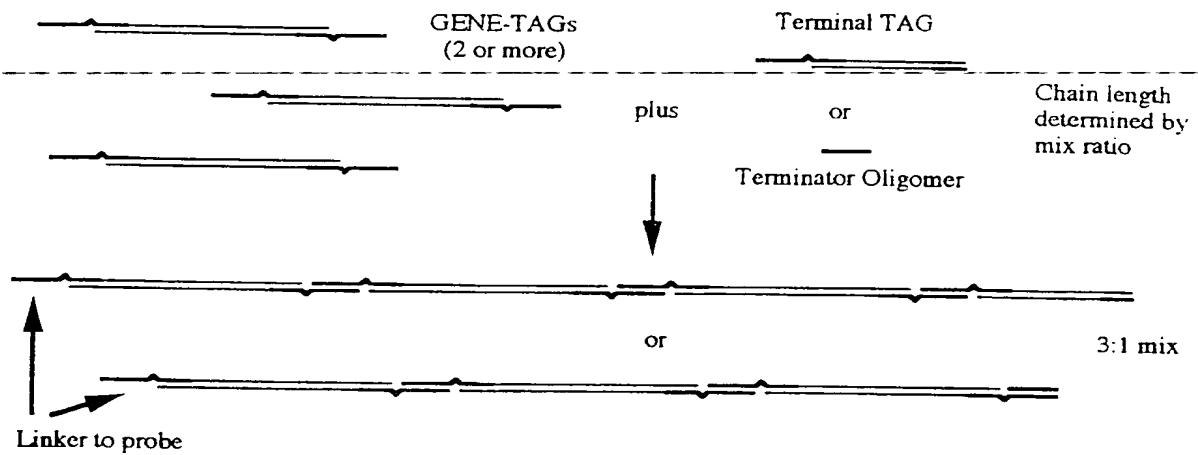B. Assemble GENE-TAG chains by ratio mix with Terminal TAGs or Terminator Oligomers



C. Hybridize GENE-TAGs to WRAP-PROBE


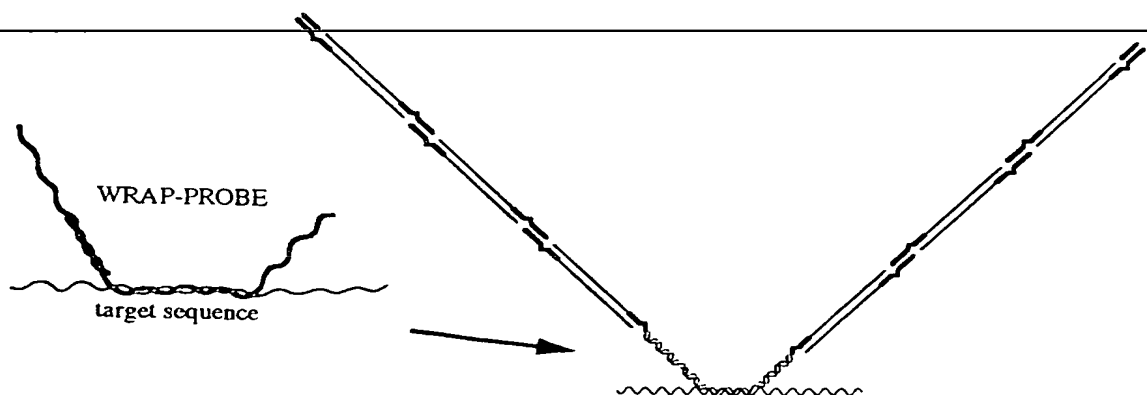


Fig. 4

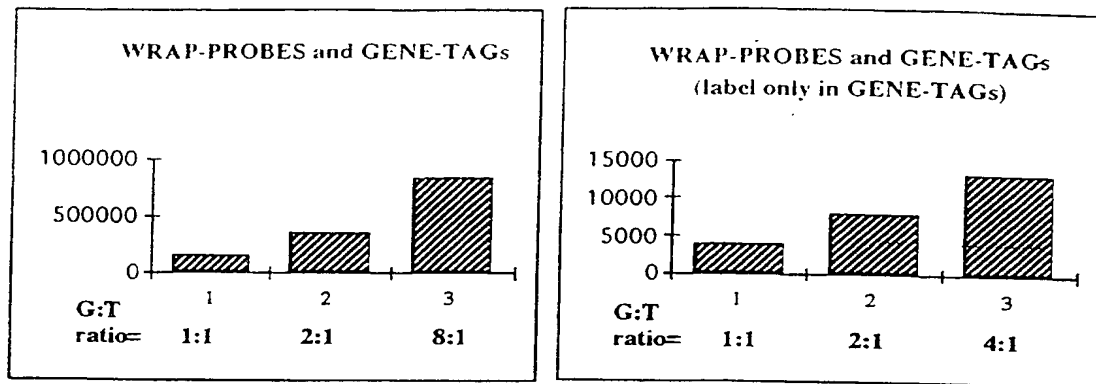A. Dot Blots of GENE-TAG Method with WRAP-PROBES (Example 4)



A1. Example: Dot blots of GENE-TAGs with WRAP-PROBE to MFB-P32-label in all TAGs. G=Gene-TAG, T=Terminal TAG

A2. Example: Dot blots of GENE-TAGs, same as Example 2A, but no label in Terminal TAGs

B. GENE-TAGs without PCR: Synthetic Oligomer Assembly


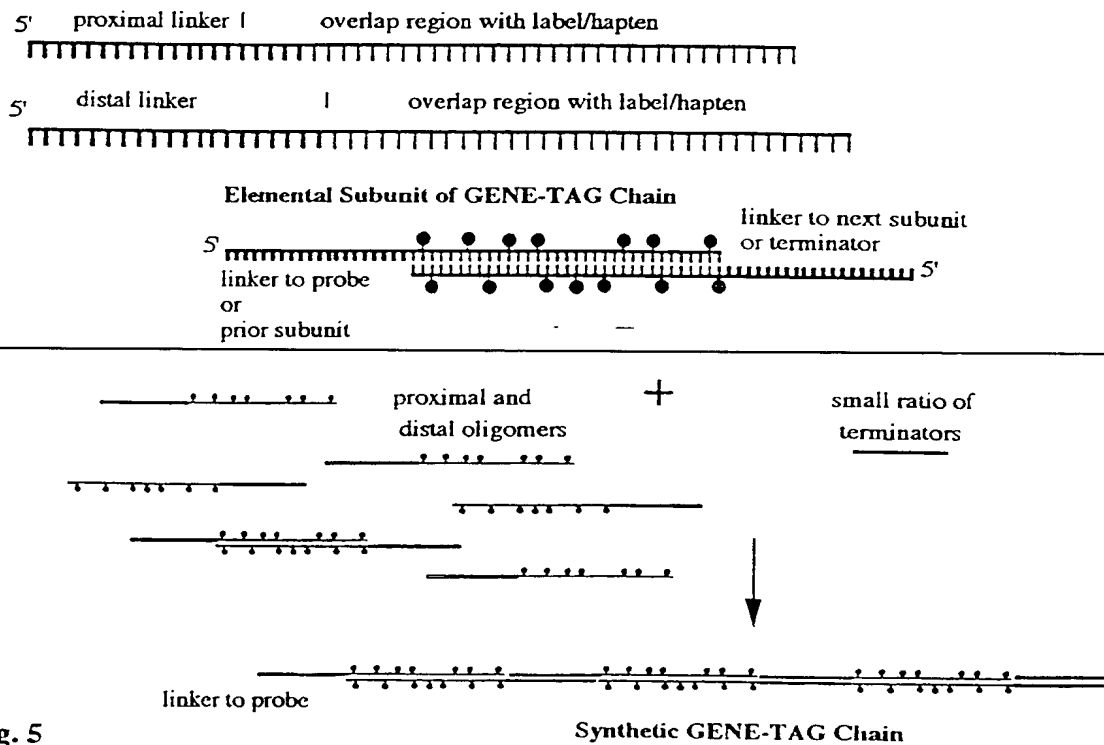


Fig. 5

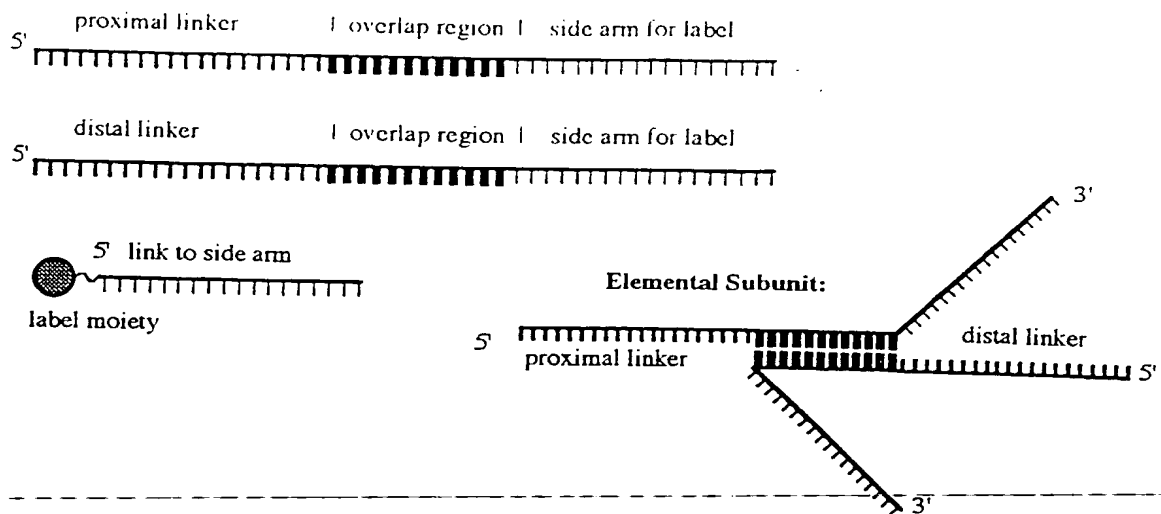A. TINKER-TAG Method: Component Design and Synthesis:



B. Component Assembly: Ratio Mix Chain Oligomers with Terminators/Label Oligomers


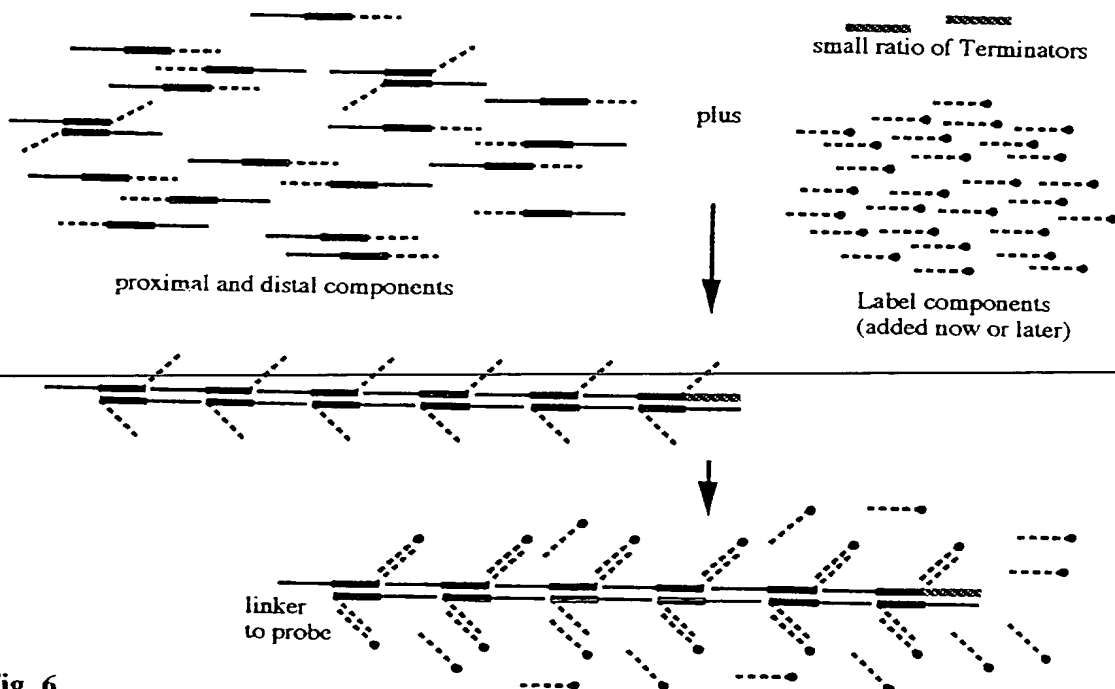


Fig. 6

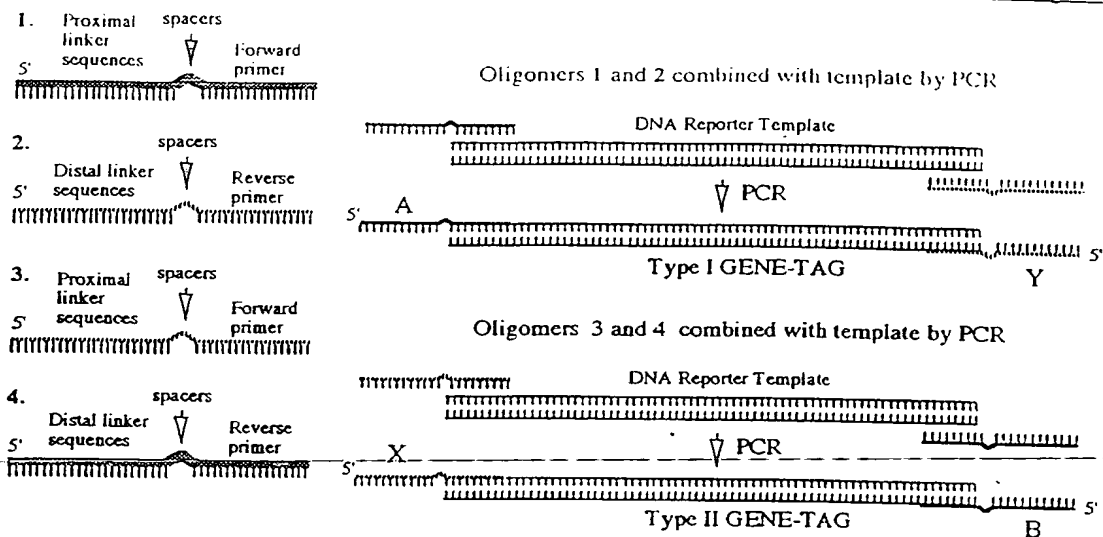A. Duo/Layered GENE-TAG Method: Construct Type I and II GENE-TAGs with alternating linkers



B. Apply Type I and Type II GENE-TAGs sequentially, washing between steps


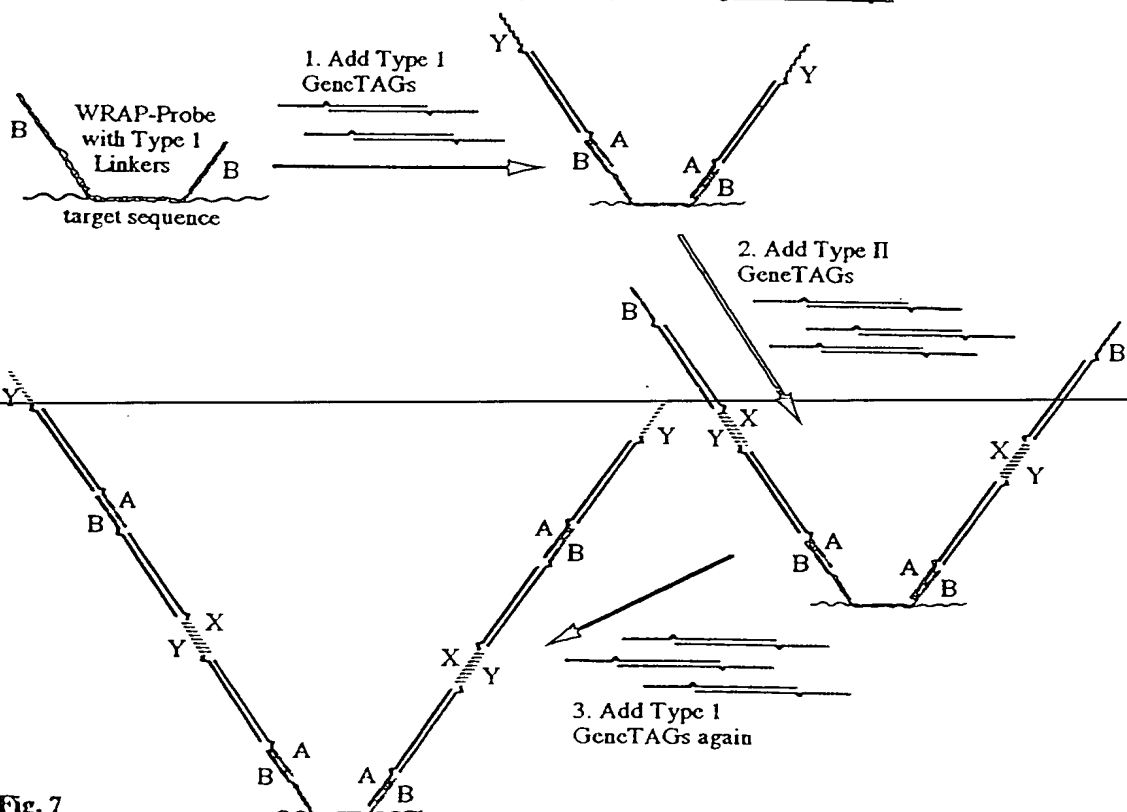


Fig. 7

A. WRAP-PROBE detection in cDNA array simulation with Chained GENE-TAGs 3:1

15-LO-1
WRAP-PROBE
applied in different
conc. then added
prejoined chain
of Gene-TAGs
of 3:1 ratio

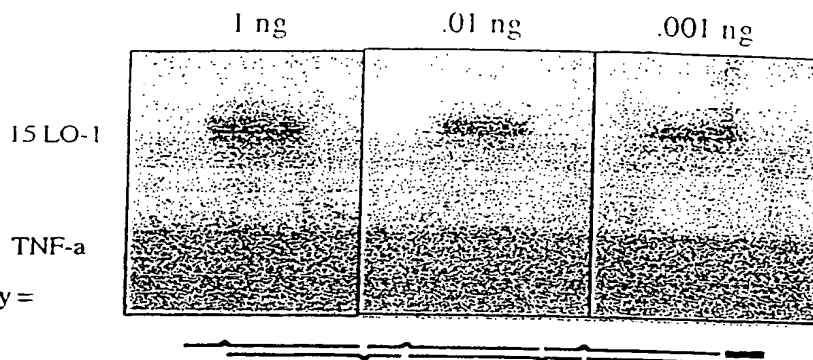Reverse dot blots on array =
 10^{11} copies per dot



B. WRAP-PROBE detection in cDNA array simulation with Layered GENE-TAGs (I+II+I)


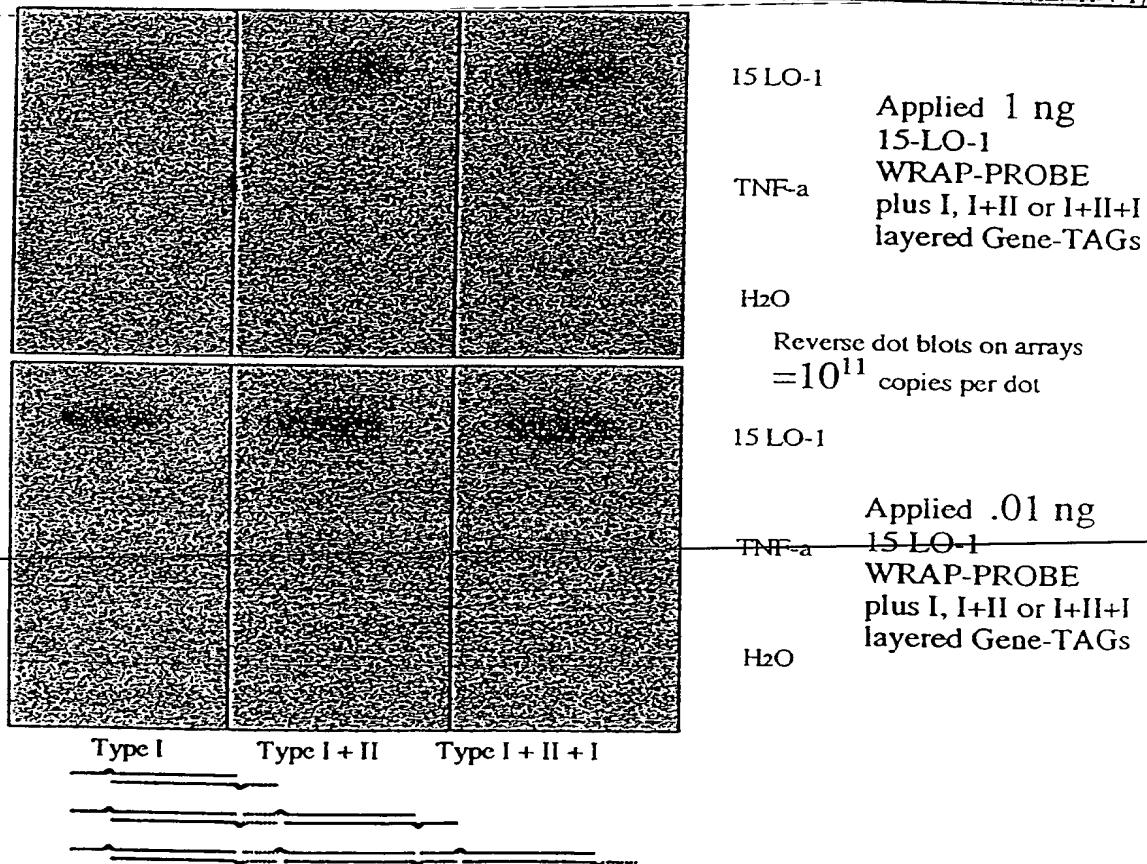


Fig. 8

Double-Duo GENE-TAG Method: Type I and Type II plus 2 distal Linkers

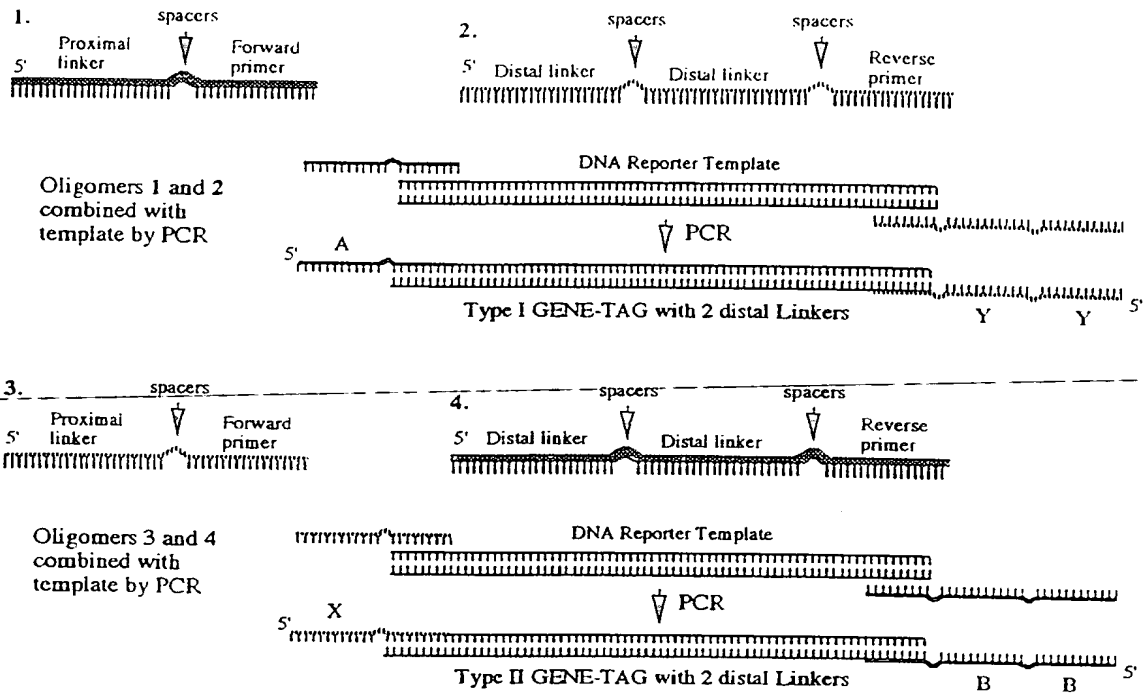A. Synthesize four GENE-TAG components and combine with reporter templates by PCR



B. Apply Double-Linker Type I and Type II GENE-TAGs in branching layers


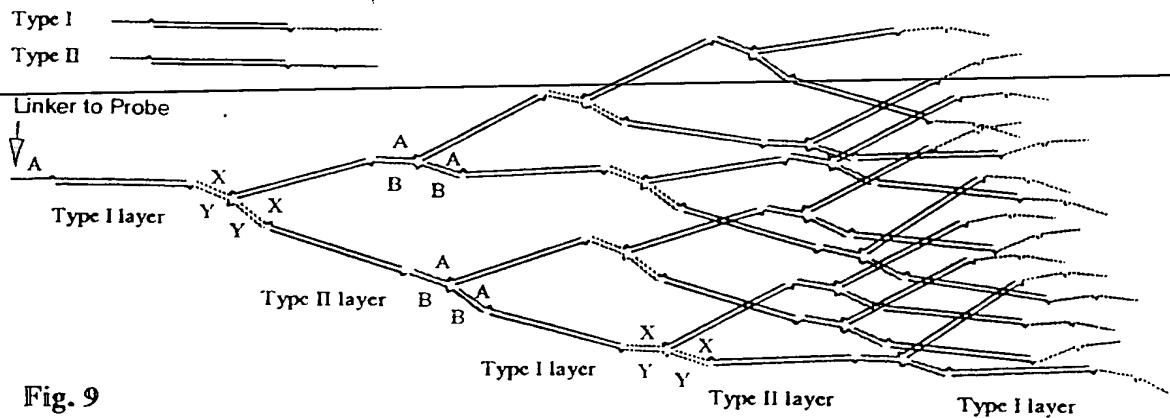
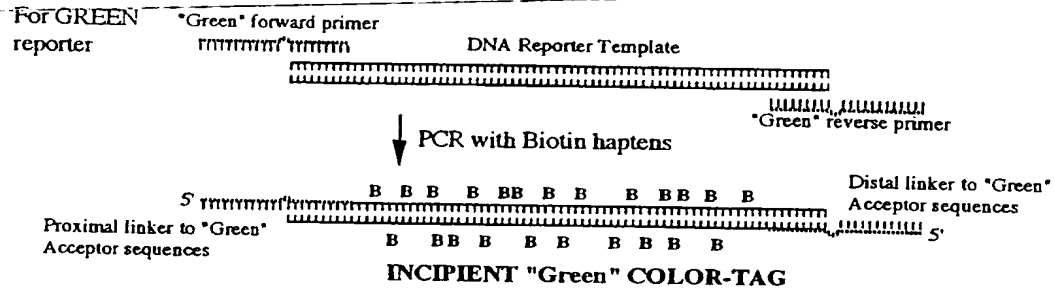


Fig. 9

COLOR-TAG Method with WRAP-PROBES based on "Red" and "Green" COLOR-TAGs

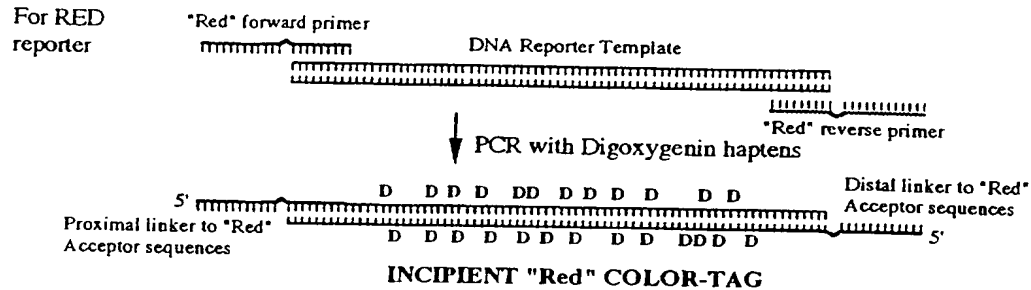A. Manufacture COLOR-TAG Reporters with different linkers and labeling



B. Construction of three WRAP-PROBES for COLOR-TAG Application:

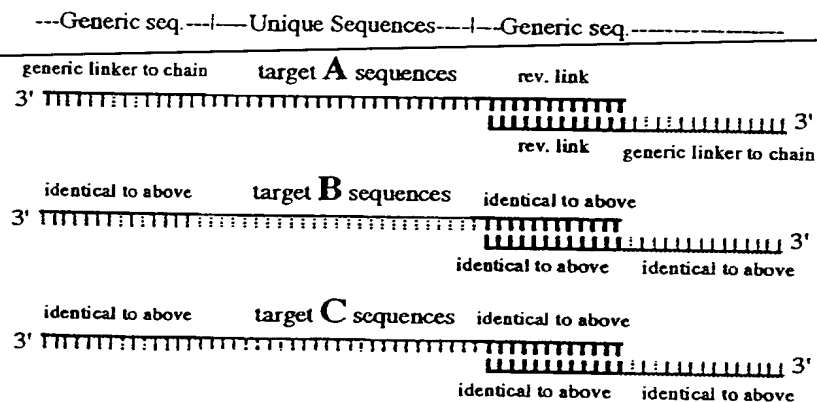


Fig. 10

COLOR-TAG Method (Part 2):

A. Synthesize COLOR-LINKER Chain Components:

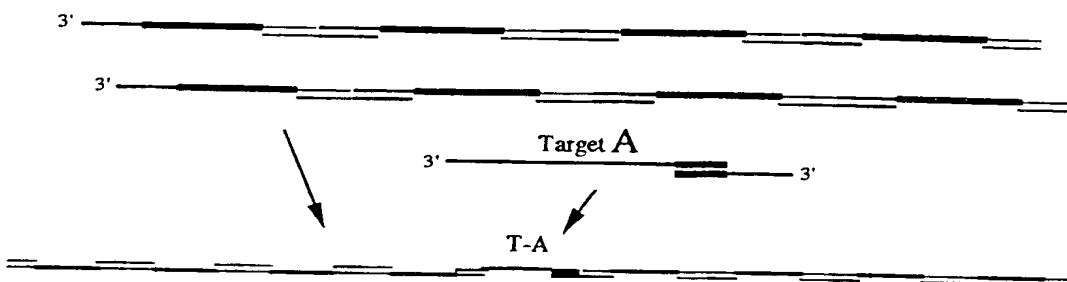
proximal R-Tag distal
linker---sequence---linker
3' ————— "Red" Acceptor Oligo 5' ————— Overlap Oligo

same as G-Tag same as
above---sequence---above
3' ————— "Green" Acceptor Oligo 5' ————— Terminator Oligo


same as B-Tag same as
above---sequence---above
3' ————— "Blue" Acceptor Oligo

Note: Acceptors have 3' end to the left.
"Blue" Acceptor is illustrated here but
is not used in example below.

B. Assemble "Red" COLOR-LINKER chain, mix on 2:1 basis with Target A Probe (eg. ABR)



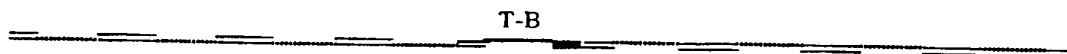C. Similarly Assemble "Green" COLOR-LINKER chain with Target B Probe (eg. D17S379)



D. Assemble 1/2 "Red" 1/2 "Green" COLOR-LINKER chain with Target C Probe (eg. CHR-12)

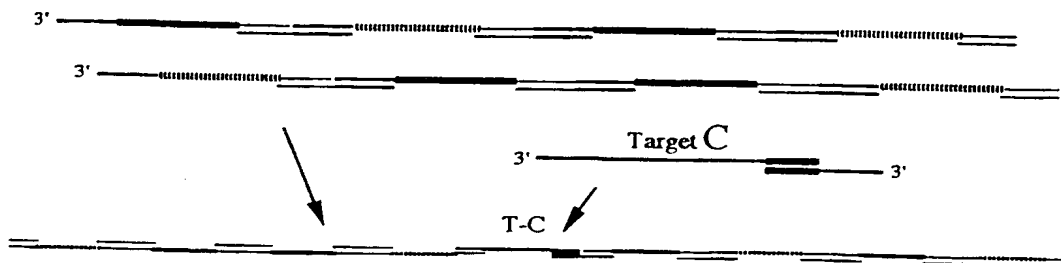
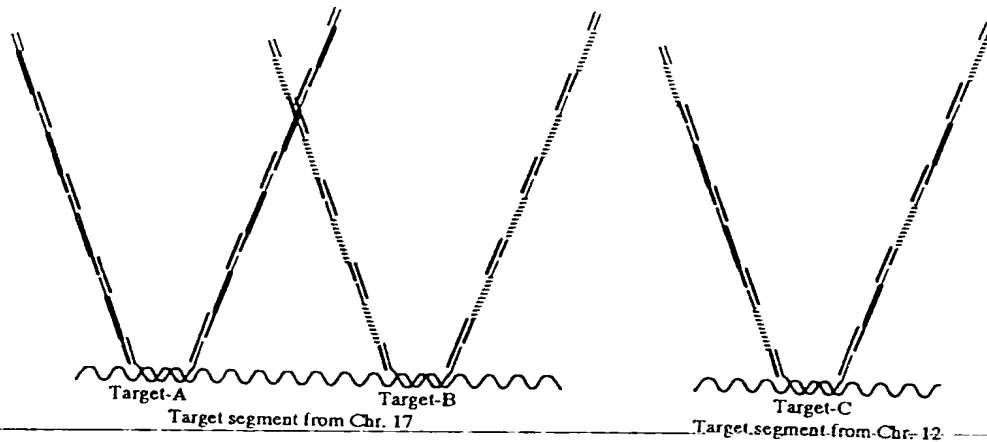
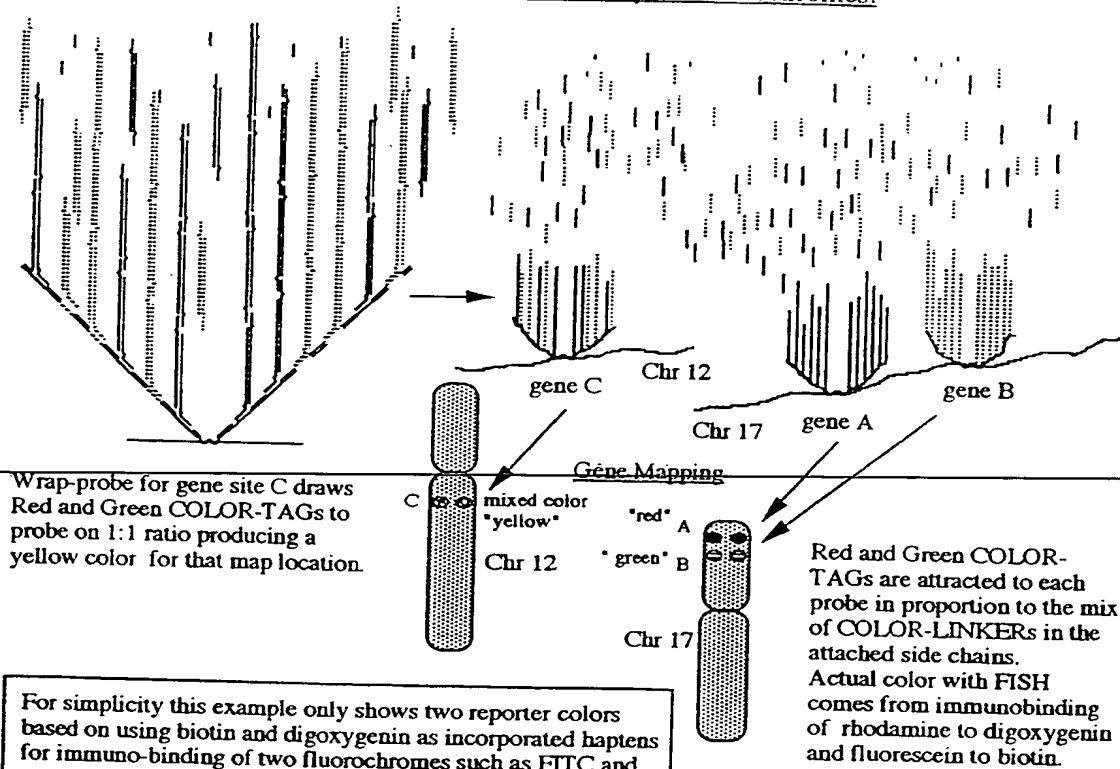


Fig. 11

COLOR-TAG Method (Part 3):A. Apply WRAP-PROBES with probe specific COLOR-LINKER chains to TargetsB. Apply set of COLOR-TAGs, plus Terminators, plus Fluorochromes:

For simplicity this example only shows two reporter colors based on using biotin and digoxigenin as incorporated haptens for immuno-binding of two fluorochromes such as FITC and rhodamine. Three color mixing can be easily achieved with alternate haptens and immuno-bound fluorochromes or with incorporation of bases that have fluorochromes directly attached.

Fig. 12

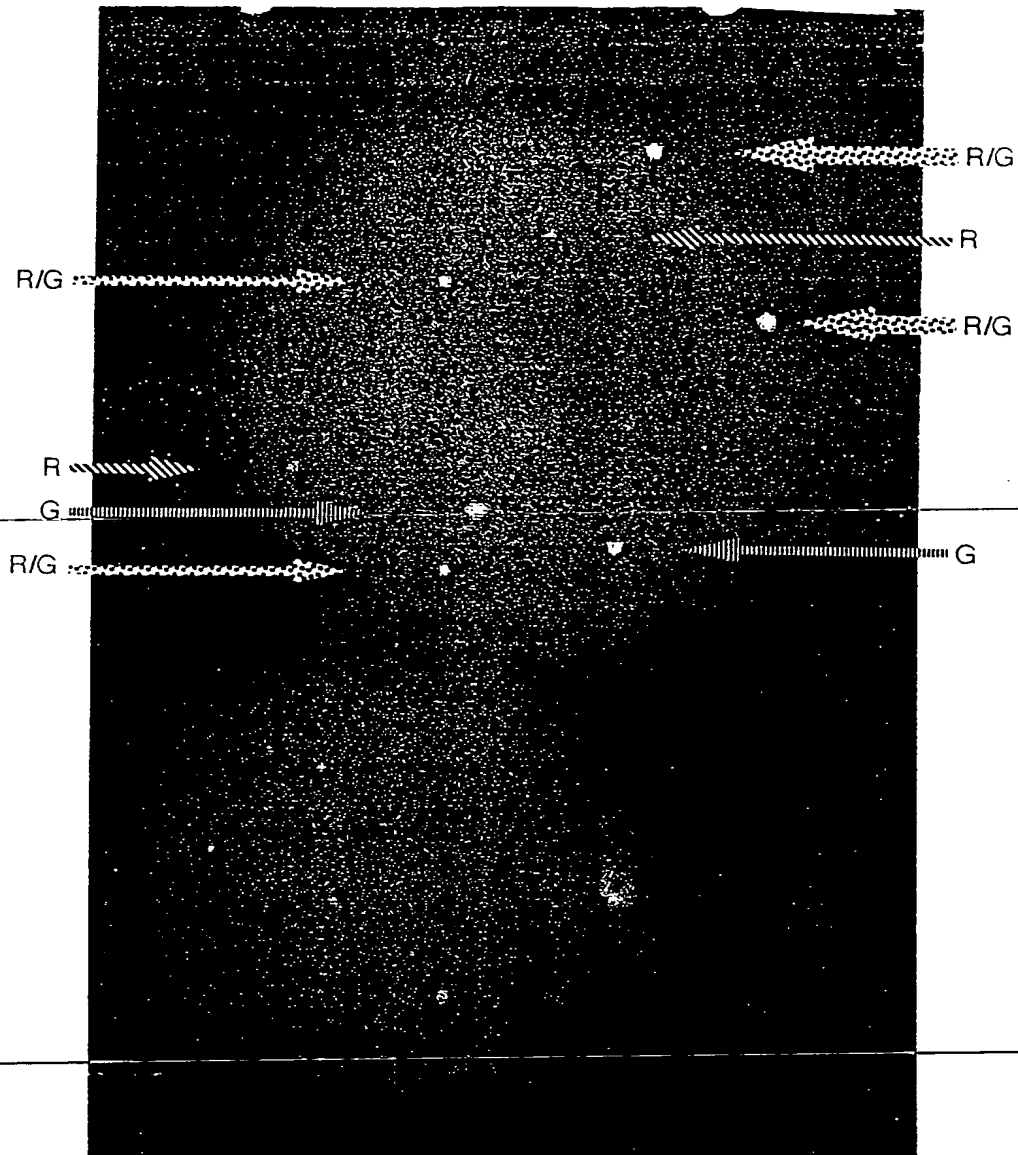




Photo image of FISH detection with four WRAP-PROBES with attached chains of either Red COLOR-TAGS, Green COLOR-TAGS or mixed Red and Green COLOR-TAGS. In nuclei, four pairs of detected dots are expected. For this gray tone copy of the color image, colors and sizes indicated by arrow patterns

Chr. 12 short repetitive site, mixed red/green R/G: 

15-LO, green: G:  ABR, red: R: 

Marker, mixed R/G:  All four WRAP-Probes have a target of 30 bp or less

Color filtering provides a clear spectral discrimination of the detection of a specific color and the presence of mixed color

Fig. 13

A. Multi-LINKER Method: Elemental form as Single Synthetic Oligonucleotide

ONE-TO-TWO Multi-LINKER

Prox. Linker 1 Dist. Linker A Dist. Linker B

Generally a Proximal probe-specific Linker (5' or 3') plus two or more Distal reporter-specific Linkers

Alternative ONE-TO-TWO Multi-LINKER with Spacers

Prox. Linker 1 Dist. Linker A Dist. Linker B

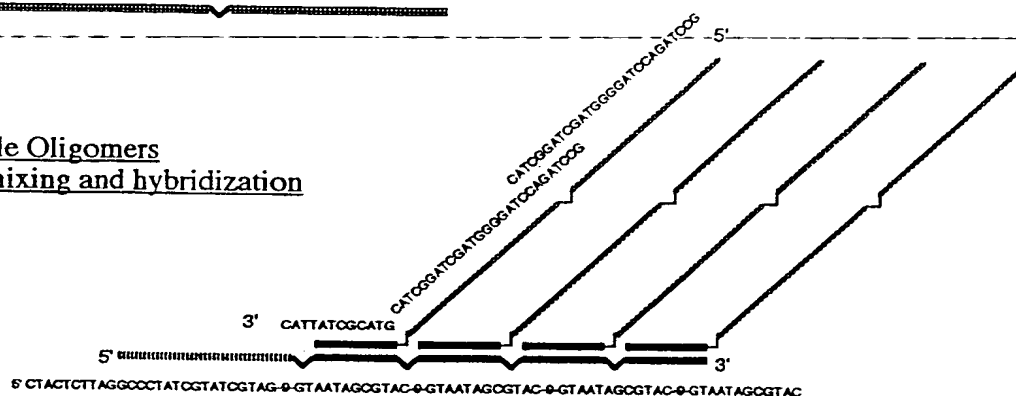
B. Preferred Two-Part Multi-LINKER embodiment for binding 8 GENE-TAGs

Synthesize Oligomers:


Five Prime ONE-TO-FOUR First Linker

Three Prime ONE-TO-TWO Second Linker

Assemble Oligomers by 4:1 mixing and hybridization



Assemble WRAP-Probe by hybridization

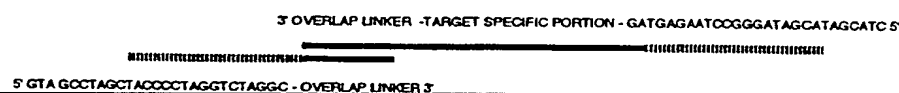

Assemble and crosslink probe and Multi-LINKER units

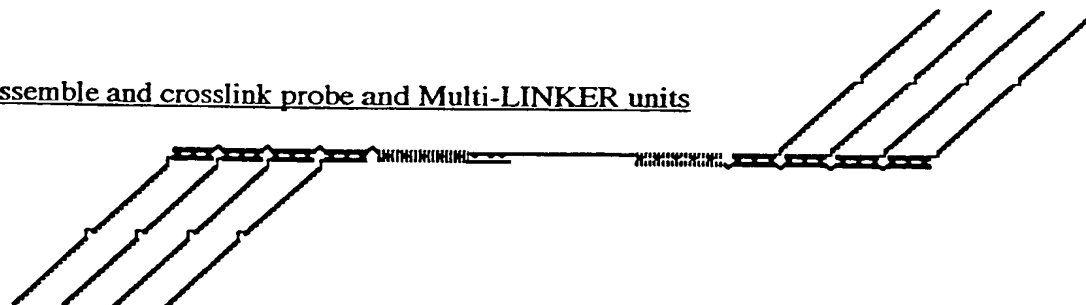


Fig. 14

A. Multi-LINKER Method: Preferred Three-Part embodiment for binding a multiplicity of short oligomers that have 5' prejoined labeling agents

Synthesize Multi-LINKER "Red" set and oligo label units with "Red" fluor (eg. Cy5)

5' CTACTCTTAGGCCCTATCGTATCGTAG-9-GTAATAGCGTAC-9-GTAATAGCGTAC-9-GTAATAGCGTAC-9-GTAATAGCGTAC

~~~~~

5' CTAGGTAGCTAG-9-CTAGGTAGCTAG-9-CTAGGTAGCTAG-9-CTAGGTAGCTAG-9-GTACGCTATTAC

~~~~~

5' CTAGCTACCTAG-99-GTACGTAAGTAG-99-GTACGTAAGTAG-99-GTACGTAAGTAG-99-GTACGTAAGTAG

~~~~~

"Red" fluor such as Cy5      5' cy5-CTAGTTACGTAC  
~~~~~

Synthesize Multi-LINKER "Green" set and oligo label units with "Green" fluor (eg. Cy3)

5' GCCTAGACCTAGGGGTAGCTAGGCTAC-9-CTACCTATCTAC-9-CTACCTATCTAC-9-CTACCTATCTAC-9-CTACCTATCTAC

~~~~~

5' CTAGGTAGCTAG-9-CTAGGTAGCTAG-9-CTAGGTAGCTAG-9-CTAGGTAGCTAG-9-GTACGCTATTAC

~~~~~ <---same as second oligo in "Red" set above

5' CTAGCTACCTAG-99-CTATCTAGTACG-99-CTATCTAGTACG-99-CTATCTAGTACG-99-CTATCTAGTACG

~~~~~

"Green" fluor such as Cy3      5' cy3-CGTACTAGATAG  
~~~~~

B. Assembly by hybridization and crosslinking

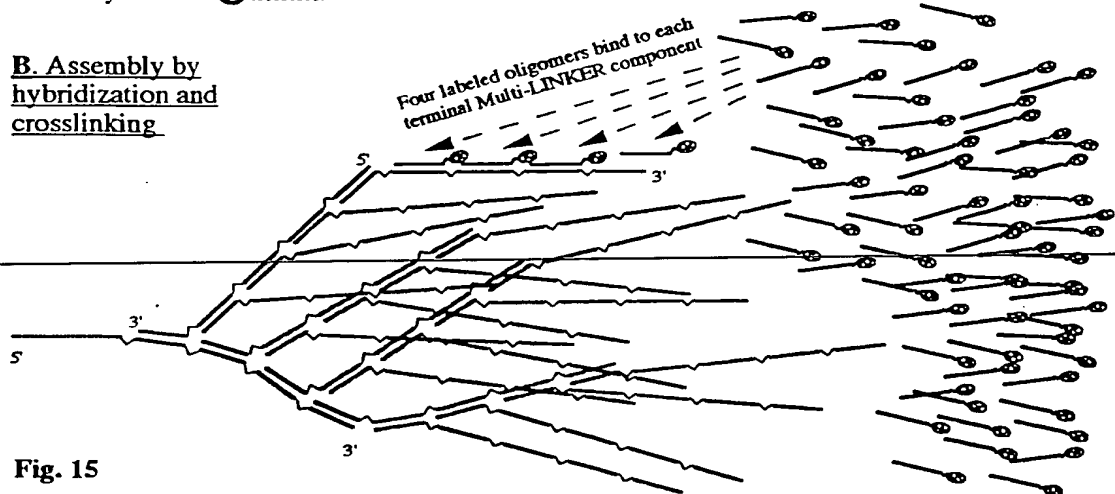


Fig. 15

GAP-LOCK Capture Probe Method: with Tinker-TAGs and Labeled Oligonucleotides

Probe components:

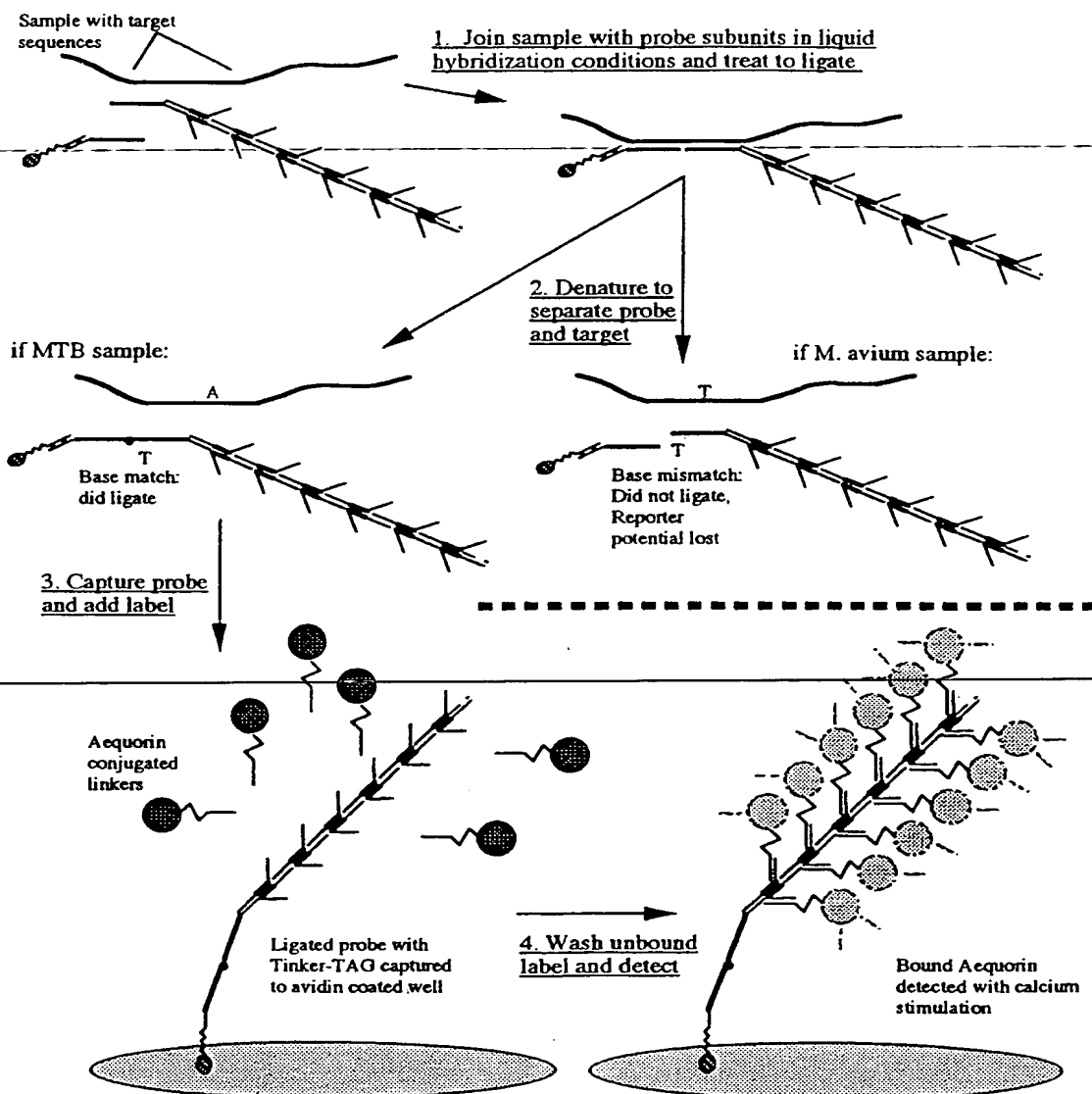
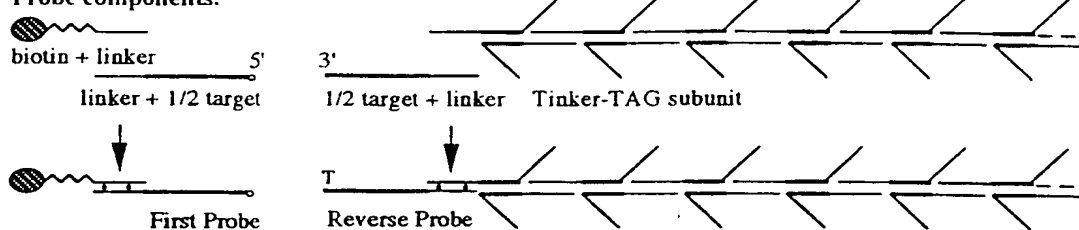


Fig. 16

RING-LOCK PROBE METHOD:A. Synthesize Target Specific GAP-LOCK Components:

- A. 5' upstream target link to 3' end of C 5' Probe subunit (First Probe)
 phosphate
- B. 5' link to 5' end of E downstream target 3' Probe subunit (Reverse Probe)

Synthesize Generic RING-TAIL Components:

- C. 5' link to reporters overlap D 3 spacers link to 3' end of A Forward Ring Subunit
- D. 5' link to reporters overlap C 2 spacers link to 3' end of E Reverse Ring Subunit
-
- E. 5' link to 5' end of B link to 3' end of D Reversing oligonucleotide

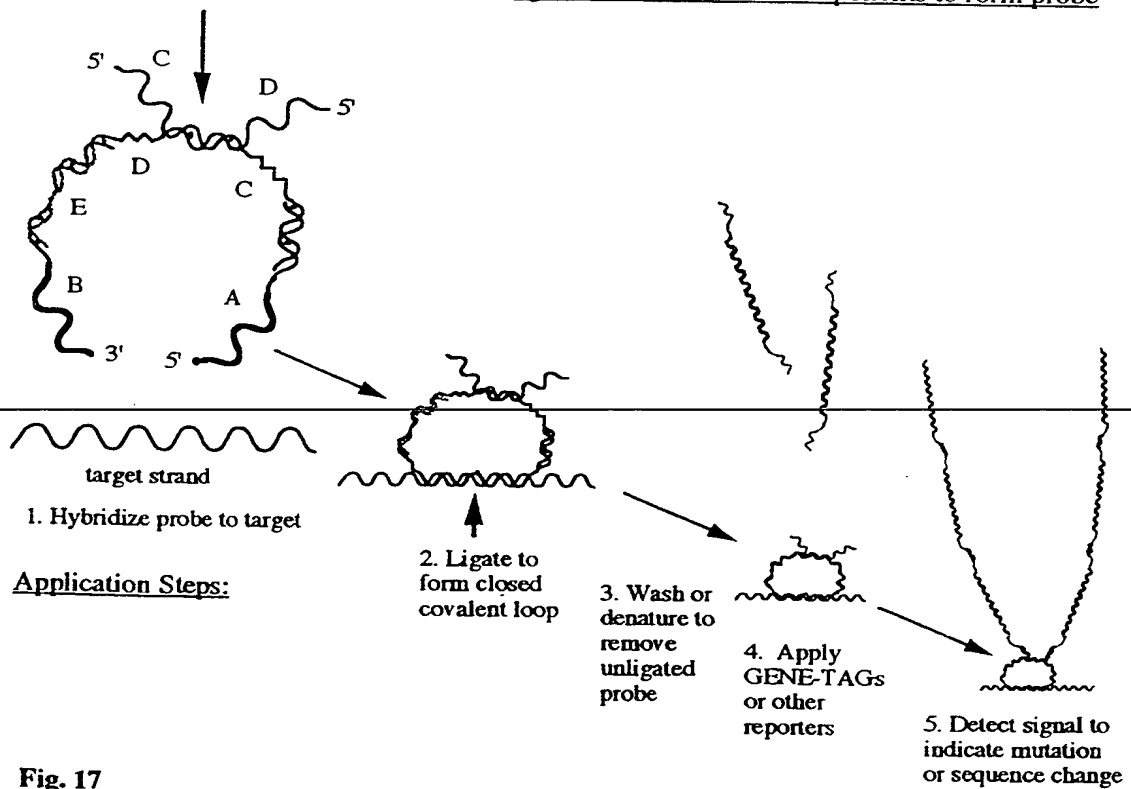
B. Assembly: Hybridize A B C D and E together and crosslink components to form probe

Fig. 17

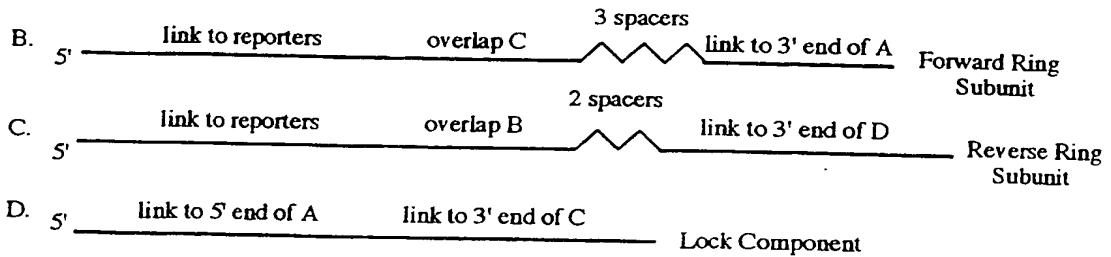
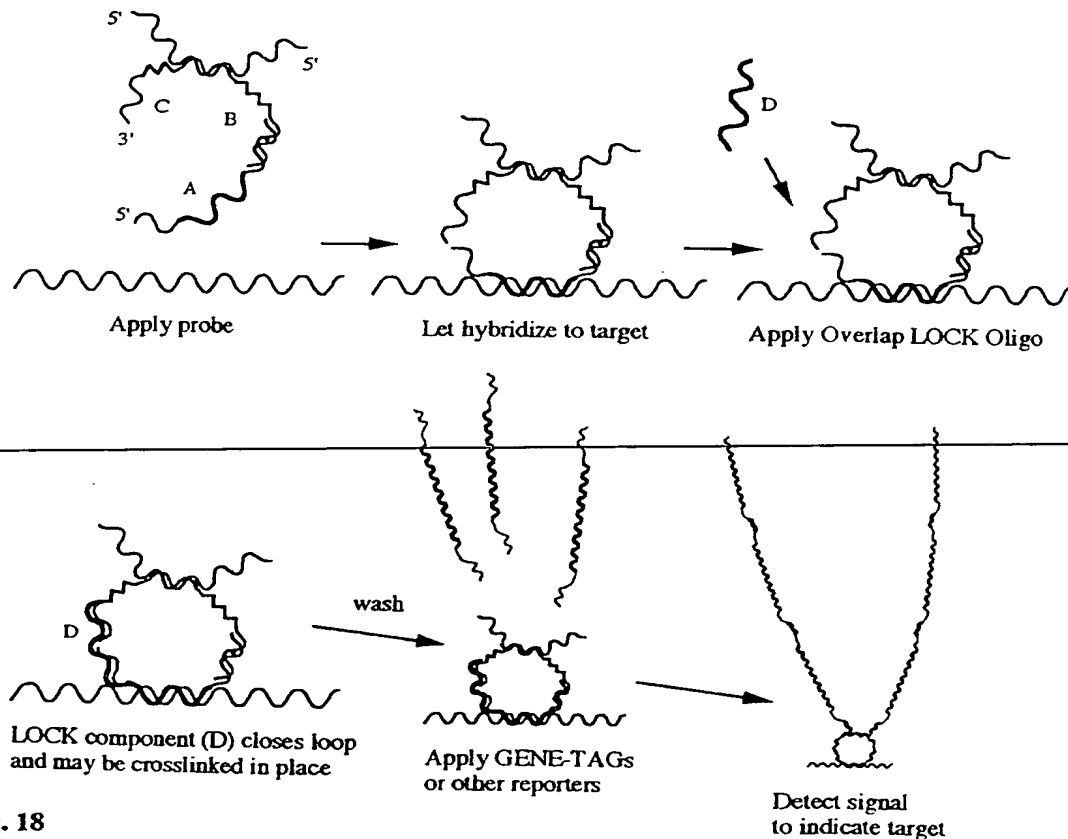
WRAP-LOCK PROBE METHOD:A. Synthesize Target Specific WRAP-PROBE Component:Synthesize Generic RING-TAIL Components:B. Assemble and crosslink four probe components and Apply to target in steps below:

Fig. 18

DOUBLE-LOCK Probe Method: Employs same probe design as RING-LOCK Probe:

- A. Components based on one probe targeting the sense strand alteration, the other targeting the complementary anti-sense strand alteration. Two hits of two color or signal type indicates confirmation of detection of rare alteration or variation in sequence. The critical base targeted is depicted as the last base on the 3' probe end of each probe.

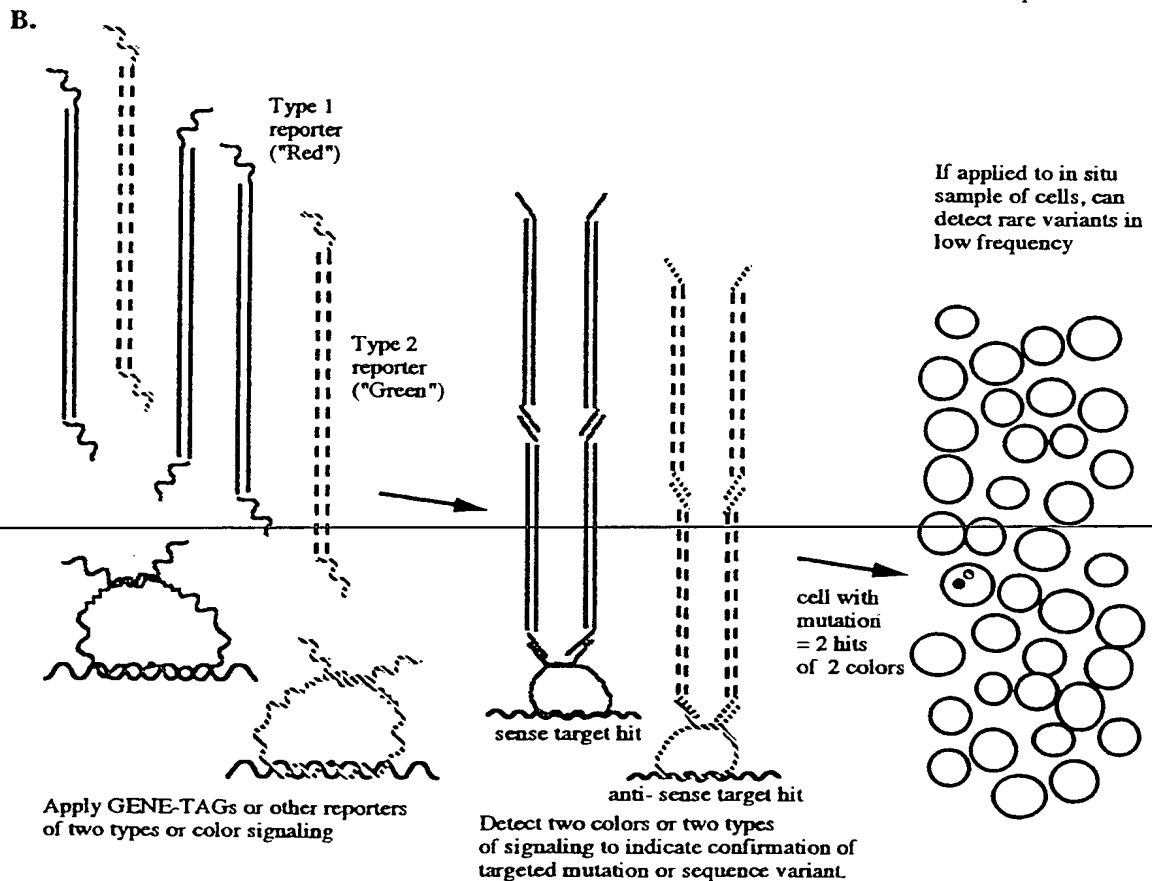
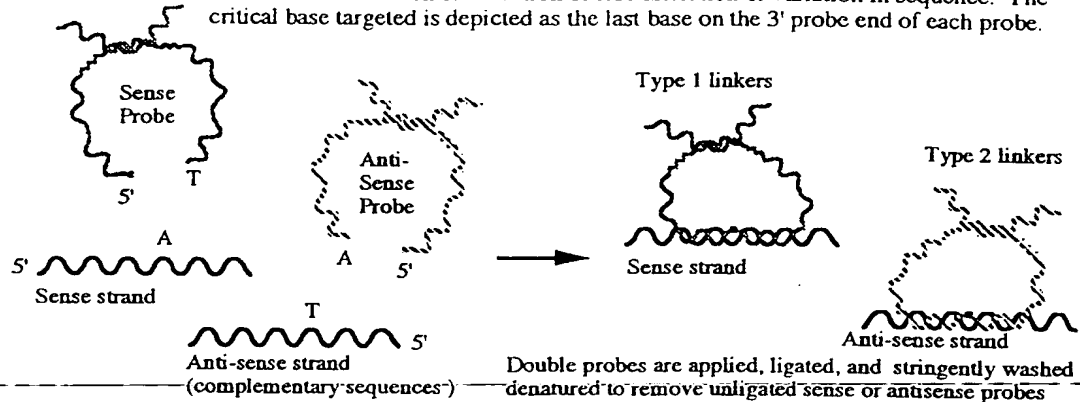


Fig. 19


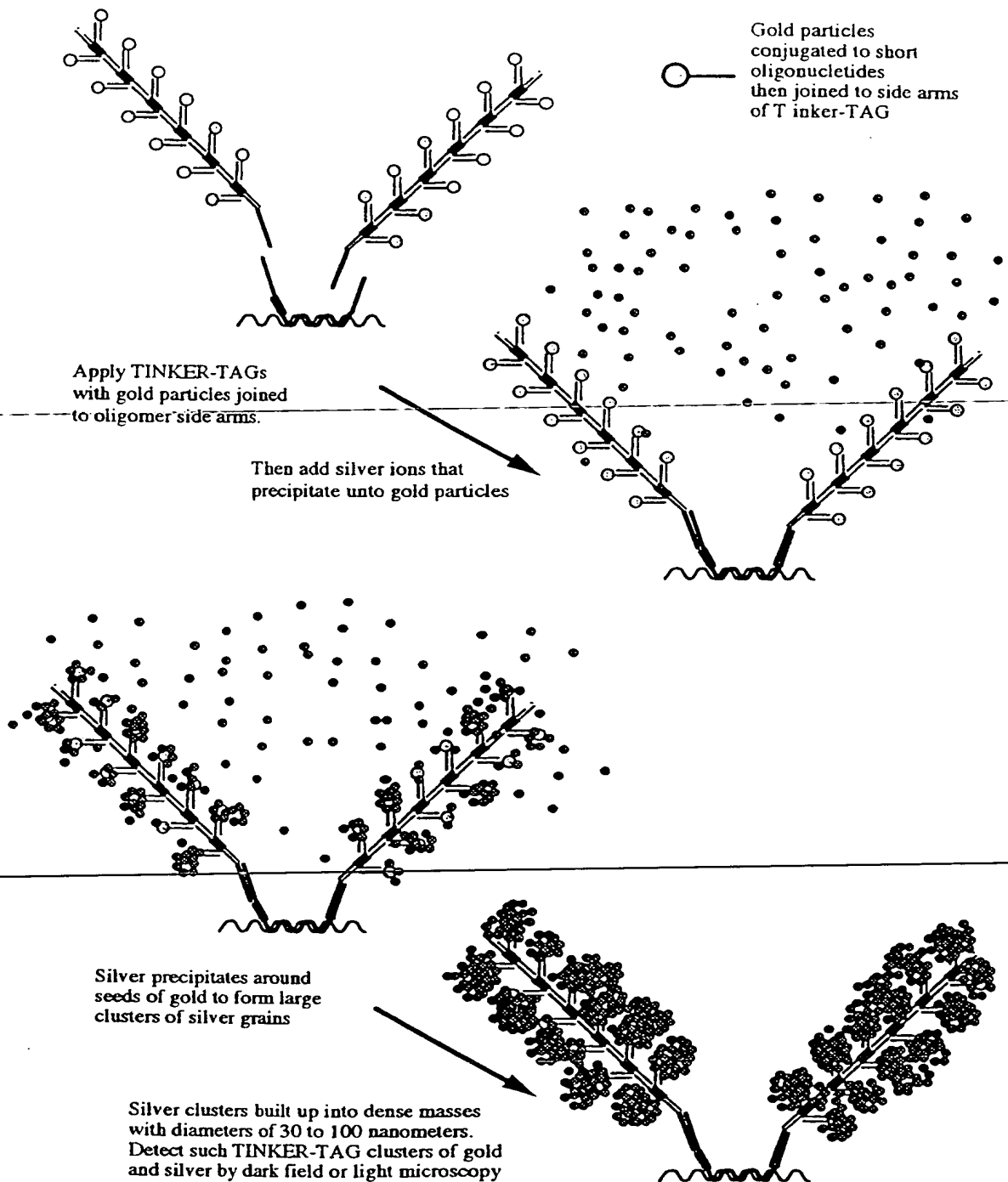
GOLD-TAG Method: Gold plus Silver TINKER-TAGs applied to WRAP-PROBES

Fig. 20

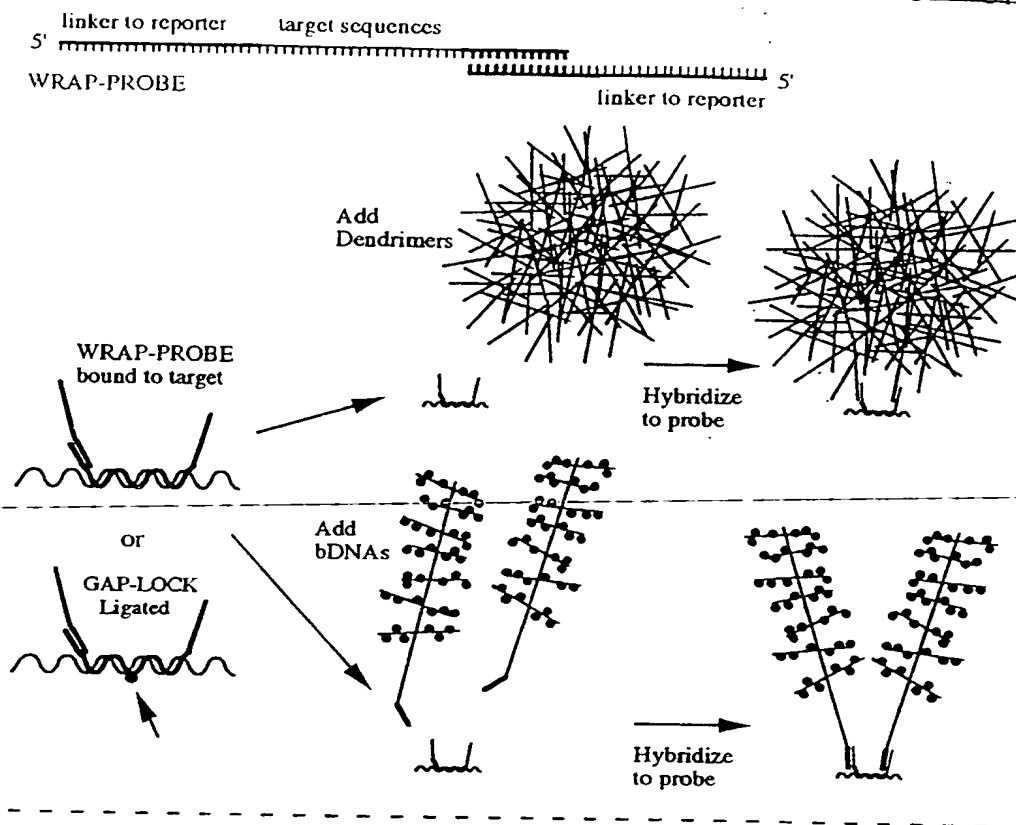A. Apply WRAP-PROBES (or GAP-LOCK) with Dendrimers and bDNA signaling products:



B. Aequorin Detection of MTB DNA using GAP-LOCK First Probe: (Example 9)

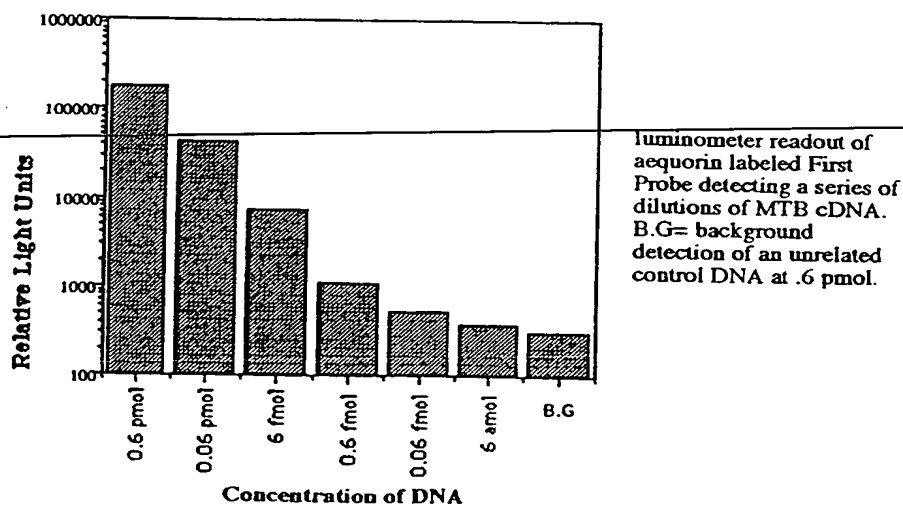


Fig. 21